

EN



ageLOC® GALVANIC BODY SPA
Owner's Manual

ageLOC® Galvanic Body Spa

Welcome to the luxury and convenience of experiencing scientifically advanced spa benefits at home. Now there's no appointment necessary for slimmer, smoother, firmer looking skin to help your body look younger.

For decades, spa and salon professionals have used gentle galvanic currents in customised treatments to refresh and energise the skin. Today, spa treatments are no longer reserved for the elite – now you can enjoy them as a part of your regular skin and body care regimen, at your leisure in the comfort of your own home.

Ergonomically designed to follow the contours of your body, ageLOC® Galvanic Body Spa features a self-adjusting galvanic current and proprietary new ageLOC® body conductive surface area designed to increase delivery of ageLOC® to your skin. The new pulsating galvanic current also aids in stimulating skin and increasing movement of fluids to purify and leave your skin feeling more refreshed than ever before.

Now you can unlock the secret to slimmer, firmer looking skin and experience the benefits of an at-home spa treatment with the ageLOC® Galvanic Body Spa.

Before you begin


Please read the following important warnings before using your ageLOC® Galvanic Body Spa instrument:

1. The ageLOC® Galvanic Body Spa instrument is only for use on healthy skin. DO NOT use over an inflamed or open wound, if you have allergies to metal, or have problematic or overly sensitive skin. Check with your physician before using the ageLOC® Galvanic Body Spa instrument if you are pregnant, have a pacemaker or similar device, are epileptic, have a metal implant, or in case of illness.
2. DO NOT use the ageLOC® Galvanic Body Spa on the face or other sensitive areas. The instrument and products are for external use on targeted areas only. If you experience any prolonged redness or abnormal irritation, discontinue use immediately. If conditions persist, consult your physician. Keep out of reach of children.
3. New ageLOC® Galvanic Body Spa users and those with sensitive skin or allergy prone skin should always test on a small area of the skin before using the product as directed.
4. Use product only for its intended use as described in the Quick-Start Guide. DO NOT use longer than the recommended times or frequencies listed in this guide.
5. Remove all protective films before beginning.
6. This instrument contains batteries. To reduce the risk of shock, burn, fire or injury: Do not immerse your instrument in water. Do not place or store where it can fall or be pulled into a bathtub, shower, sink or toilet. Never operate the instrument if it has been damaged. Do not use the instrument in any way not described in the Quick-Start Guide.
7. Please read through this manual in its entirety prior to using your ageLOC® Galvanic Body Spa instrument or accompanying products.

To ensure the correct waste treatment of the ageLOC® Galvanic Body Spa instrument and battery, please dispose of them in accordance with any applicable local laws or requirements for disposal of electrical equipment/batteries. Should not be disposed of as general household waste.



Your ageLOC® Galvanic Body Spa system includes the following:

- 1 Carrying Pouch
- 2 Quick-Start Guide
- 3 CD-ROM with full Owner's Manual
- 4 ageLOC® Galvanic Body Spa instrument (with two AAA batteries installed)
- 5 Instrument Stand

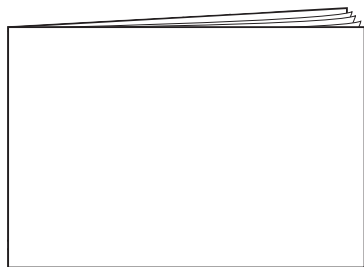

2




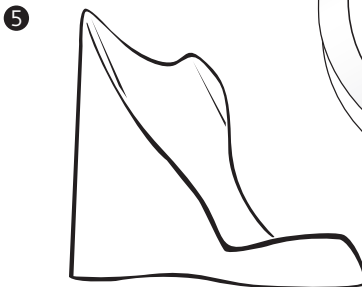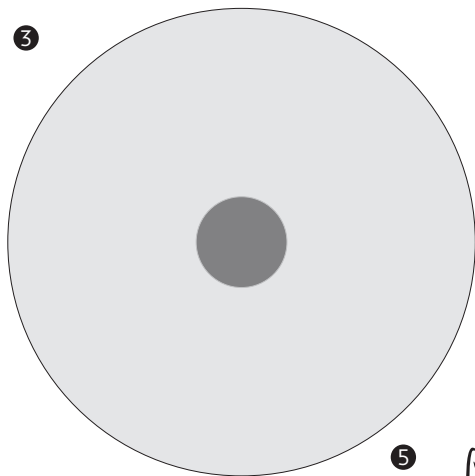

Note: The products designed for use with the ageLOC® Galvanic Body Spa are sold separately.

Your ageLOC® Galvanic Body Spa instrument has the following features:


- 1 Display
- 2 On/Off Button
- 3 Chrome Activators
- 4 Battery Compartment (two AAA batteries already installed)
- 5 ageLOC® Conductive Surface

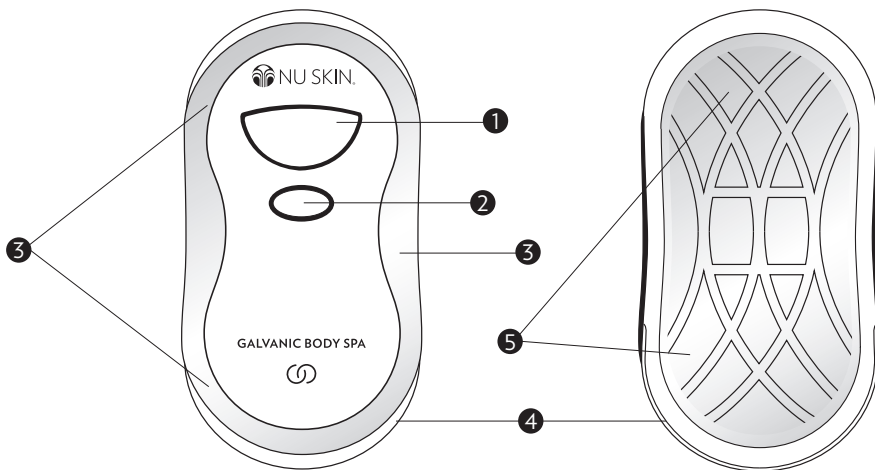

Your display has the following icons:

- 1 Low Battery Indicator – this icon will illuminate when batteries are low.
- 2 ageLOC® rings – this icon appears while the unit is on.
- 3 Clock face and number 5 – these icons will appear to indicate the length of treatment time, which is 5 minutes.

Note: Remove protective plastic film from instrument before use.

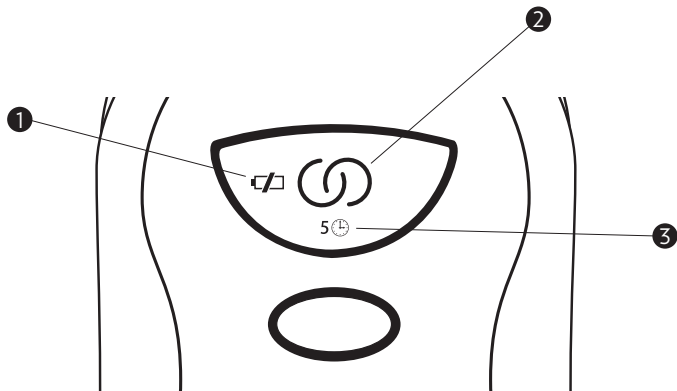

How to use your ageLOC® Galvanic Body Spa instrument

Backlight: The backlight will remain lit while the instrument is on. If the backlight is blinking, it is either because the instrument has lost contact with your skin or there is not enough moisture for optimal conductivity. To continue your spa treatment, ensure adequate moisture at all contact points and resume contact with your skin.



Product application: Apply the appropriate gel or treatment formula liberally to the area being treated.

Treatment selection: Press the on/off button once to turn on the instrument. The ageLOC® rings appear in the middle of the display and the treatment time of 5 minutes appears below the ageLOC® rings. In order for a treatment to begin, gently touch the conductor's surface against the treatment area for the current level to self-adjust. The instrument will beep one, two, or three times, indicating that it has self-adjusted to your skin.

To activate the galvanic current: Your moistened hand must be touching the chrome activators while the conductor makes contact with the treatment area. Make sure that all contact points are thoroughly moistened. Remember, all contact points must be wet to optimally conduct the galvanic current.

During treatment: The instrument will beep every 10 seconds to indicate the unit is on and working. Once treatment time reaches 30 seconds remaining, it will beep every 5 seconds until treatment is complete.

Treatment interruption: To maintain the galvanic current throughout the treatment, you must keep the conductive surface in contact with the area being treated and ensure adequate moisture at all contact points. Breaking contact with the treatment area or chrome activators will disrupt the current and interrupt the treatment. In this case, the backlight will flash on and off and the instrument will stop beeping until the circuit is either re-established or until the instrument shuts off after two minutes of non-use.

Treatment completion: At the end of your treatment, the clock and number of minutes will disappear, but the ageLOC® rings will remain on the screen. A “treatment complete” sound will also play simultaneously. If you wish to initiate another galvanic treatment, press the on/off button and re-initiate the cycle. To turn the instrument off, simply press the on/off button twice.



ageLOC® Body Shaping Gel

For optimal results, use ageLOC® Body Shaping Gel with your ageLOC® Galvanic Body Spa two to three days each week, morning or night, on legs, arms, and/or lower abdominal area. Read through these instructions completely before starting your body shaping treatment.



1. Apply ageLOC® Body Shaping Gel:

- Smooth gel liberally on area to be treated. If you plan to treat more than one area, only apply gel to the area to be immediately treated.
- Press the on/off button once. The ageLOC® rings will appear in the middle of the display. Treatment time of 5 minutes will also appear below the ageLOC® rings.
- Check to ensure that your fingers or hand are adequately moistened and in contact with the chrome activator areas.
- Touch the conductor head against your skin for the current level to self-adjust.

2. Perform a treatment with ageLOC® Body Shaping Gel:

First, thoroughly moisten your fingers and the palm of your hand with water or NaPCA Moisture Mist. Hold the instrument with an adequately moistened hand, making sure your fingers are touching the chrome activator. Ensure that your fingers do not touch the conductive surface on the underside of the instrument, as this will not allow the instrument to self-adjust to the skin on the targeted area. Touch the conductive surface to your skin where you have applied ageLOC® Body Shaping Gel. You will hear one, two or three beeps as the instrument automatically adjusts the current level for your skin.

A. FOR USE ON THE LEGS


Begin lightly massaging in a circular motion where the product has been applied, moving up the thigh and gliding back down.

B. FOR USE ON THE UPPER ARM

Begin lightly massaging in a circular motion where the product has been applied, moving up the upper arm and gliding back down.

C. FOR USE ON THE LOWER ABDOMINAL AREA

Begin lightly massaging in a circular motion where the product has been applied, moving up across the lower abdominal area.

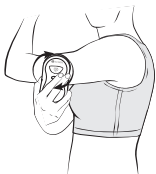

ageLOC® Body Shaping Gel (continued)

- 3. If you lose contact with the chrome activators** at any time during your treatment, simply adjust your hold on the instrument and resume your treatment. The instrument will beep every 10 seconds to indicate that it is on and working properly. Once treatment time reaches 30 seconds remaining it will beep every 5 seconds until treatment is complete.
- 4. Completing a treatment**
 - a. At the end of your treatment, the instrument will emit a long beep followed by two short beeps indicating that treatment is complete. The clock and treatment time will disappear and the ageLOC® rings will stay on for 30 seconds of inactivity. Remove excess product from your skin. Clean remaining product from the instrument and conductive surface with a damp towel and dry thoroughly.
 - b. To begin another galvanic treatment, press the on/off button and re-initiate the cycle. Follow the directions above for the area you wish to treat.
- 5. For an extension of benefits**, use ageLOC® Dermatic Effects each day, morning and night. Read your ageLOC® Dermatic Effects directions completely before using.

Note: 5 minutes is the recommended treatment time for one area, but you may not need a full 5 minutes.

Important maintenance and warranty information

Directions for changing batteries: Your ageLOC® Galvanic Body Spa comes equipped with two AAA batteries. When the batteries are low, the Low Battery indicator will appear on the display screen. When this happens, install two new AAA batteries.

To replace batteries, apply pressure just above the battery compartment cover. Grip the cover tightly and pull firmly. If the cover is difficult to remove, apply more pressure to the instrument and pull on the cover with greater force. The cover will eventually disengage from the compartment. Upon removing the battery compartment cover, remove the water tight stopper by pulling outward on the tab.

Directions for cleaning your instrument: After treatment, use a damp cloth to remove any excess product from the instrument. Do not use hard chemicals cleaning products or solvents as they may damage the product. Dry thoroughly. Store in a dry place. Do not expose to extreme temperatures or high humidity.

To ensure the correct waste treatment of the ageLOC® Galvanic Body Spa instrument and battery, please dispose of them in accordance with any applicable local laws or requirements for disposal of electrical equipment/batteries. Should not be disposed of as general household waste.

Limited two-year warranty: Nu Skin® guarantees the ageLOC® Galvanic Body Spa to be free from defects in materials and workmanship for a period of two years from the original date of consumer purchase. This warranty does not cover damage to the product as a result of misuse or accident. If the product becomes defective within the two-year warranty period, please call your local Nu Skin® office to arrange for a replacement.





Imp./Importálja/ Импортёр/ Импортёр:

Nu Skin Netherlands B.V. Marinus Dammeweg 61, 5928 PW Venlo, The Netherlands

HU: Nu Skin Eastern Europe Kft., Alkotás utca 48-50, 1123 Budapest.

RO: Nu Skin Enterprises S.R.L., Str. Cobălcescu 46, 010196 Bucureşti.

RU: 125047 Россия, Москва, 4-ый Лесной пер. д.4. Тел.: 8-800-700-19-84

TR: Nu Skin Cilt Bakımı Tic. Ltd. Şti. Ayazağa Mah. Meydan Sok. No:28 Beybi Giz Plaza Kat: 26, 34396 Maslak, İstanbul tarafından ithal edilmiştir.

©2012 NSE Products, Inc. Provo, Utah, 84601. USA • www.nuskin.com



CE



www.nuskin.com